


PLATELET ROTATOR

Robust durable design for your economic choice storage frame 360° smoothly and steadily. Platelet bags can be loaded easily on the frame basket housing. Complete air circulation across both surfaces the bag.

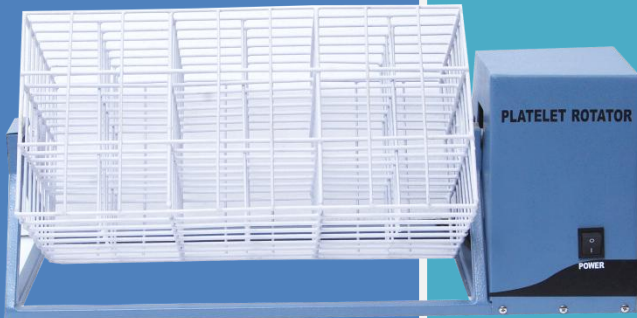

SPECIFICATION:

Model No. PR 950

Rotation Speed 6 RPM

Maximum capacity 32 Bags

Operating Volt 230V AC

Power Frequency 50/60 Hz

Power Consumption 15 watts

PLATELET AGITATOR

Platelet Agitator specially designed to maintain platelet in continues motion at a fixed speed of 70 RPM. Gentle side to side motion (1.5" or 38 mm) with 65 ± 5 strokes per minute. Removable drawers for storage of any size bags. Sturdy one piece drawers with holes for complete air circulation across both surfaces of platelet bags. Heavy duty ball bearing gear motor for noiseless and continuous operation for 24 hours a day 365 days a year.



SPECIFICATION:

Model No. PA 920

Rotation Speed 70 RPM

Maximum capacity 32 Bags

Operating Volt 230V AC

Power Frequency 50/60 Hz

Power Consumption 45 watts

Hi-Tech Instruments
Office 7 & 8nd Floor Block C-1,
Taj Medical Complex
M.A. Jinnah Road
Karachi Pakistan.
Phone: 0092 21 2787496
Mobile: 03002343299
E-Mail: hitech.aslam@gmail.com

